

Gerbera Sequencing kit v3.1

Direct substitute for BigDye® Terminator kit v3.1

Compatible with all ABI Genetic Analyzers

- No changes to workflow, protocols, calibration or settings
- Developed for use with Dye Set Z

Uniform peak heights with extremely low background

- Long sequence reads with high confidence scores
- Results that are identical to BigDye® Terminator v3.1


Gerbera Sequencing kit v3.1 contents:

- 2.5x cycle sequencing mix
- 5x Sequencing Buffer
- 125, 250, 500 and 1000 reaction packages



Can you tell the difference?

Upper panel – Gerbera Sequencing kit v3.1, lower panel – BigDye® Terminator kit v3.1

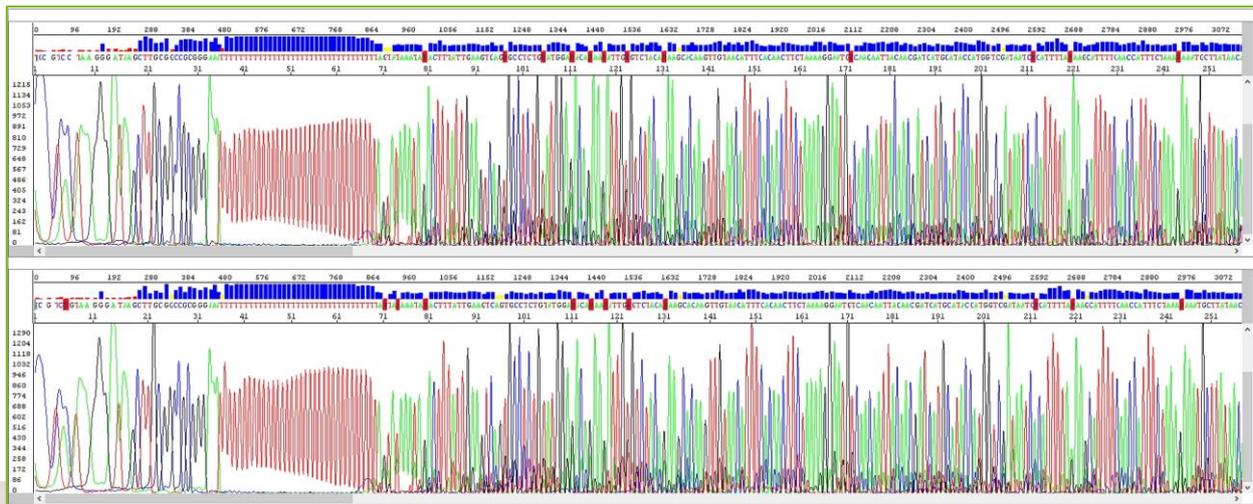

Gerbera-7 Polymer

Direct substitute for POP-7™ polymer

Preformulated separation matrix for capillary electrophoresis

- Optimized for short- to long-read sequencing and fragment analysis
- Easy to start with – you can use the current spectral calibrations

Reproducible Results

- Gerbera-7 Polymer dynamically coats the capillary wall to control electro-osmotic flow.
- Its quality and uniform consistency help to ensure reproducibility.



Packaging 7 or 28 ml

Custom packaging
available on request

